

Specifications

System Function	Category	Random access, fully automatic chemistry analyzer	
	Throughput	240 tests/hour(without ISE),315 tests/hour(with ISE)	
	Methodology	End point, Fixed Time, Kinetic	
	Reagent	Open system	
	Reagent Position	40	
	Sample Position	40	
	Minimum Reaction Volume	150uL	
	Carryover	<0.005%	
	Water Consumption	Low to 5L/Hour	
	Maximum Reaction Time	19.1 minutes (Single reagent), 15.7 minutes(Dual reagent)	
	Maximum Reaction Volume	500 uL	
	Principle	Colorimetry, Turbidimetry	
	Testing Mode	Regular mode(single&dual reagents), fast mode(single reagent)	
	Sample Mode	Random access, STAT sample priority	
Sample Reagent Probe & Mixer	Sample Volume	2~50uL, step by 0.1 uL	
	Sample Disk	40 positions	
	Sample/reagent Probe	One probe for sample and reagent with liquid level detection, vertical and horizontal collision protection and reagent volume real time monitoring function	
	Sample Tube	Regular sample tube (13 x 100mm, 13x75 mm, 12x100 mm, 12x75 mm), Vacuum tube, eppendorf tube, plastic tube, etc.	
	Reagent Volume	R1:150-450 uL, R2:10-300uL, step by 1uL	
	Reagent Disk	40 reagent positions, 24 hours water cooling, 4-12°C	
	Reagent Bottles	Compatible with Olympus and Hitachi reagent bottles	
	Mixer	High polished nano material mixer	
	Reaction System	Reaction Cuvettes	81 cuvettes, semi-permanent rigid UV special plastic cuvettes
		Wash Station	6-step auto wash station
		Reaction Temperature	37±0.1°C
		Heater	Metal thermostat for reaction disk
	Optical System	Lamp	Halogen Tungsten lamp, 12V20W, about 2000hours, water cooling system
		Filter	FMSS(Full-sealed Matrix Spectrometric System)
Wavelength		8 wavelengths: 340,405,450,510,546,578,630,670nm	
Resolution		0.0001Abs	
Linearity Range		0-3.5Abs	
Accuracy		0.5A: < ± 0.02AbS, 1 .0A: < ± 0.04Abs	
Stray Light		>4.5	
Stability		<0.01Abs/hour	
CV		<1.5%	
Half Band Accuracy		≤ ±2nm	
Cal&QC		Detector	Photodiode detector array
		Calibration Mode	One point linear, two point linear, multi-point linear, Logit-4P, Logit-5P, spline, exponential, polynomial
		Calibration Curve	Calibration curve auto check, auto curve fitting
		QC Rule	Westgard multi-rule. Cumulative sum check, Twin Plot
Software	QC Curve	Westgard multi-rule, Cumulative sum check, Twin Plot	
	QC Warning	Out of control auto warning, data automatic record and analysis	
	Operating System	Windows 7/10	
	Data Storage	Decided by computer hard disc capacity	
	LIS Interface	Bi-LIS interface	
	Print Mode	Multi default formats, user-defined formats	
	System Monitor	Real time monitoring for sample disk, reagent disk, reaction disk. Real time monitoring for QC status. Real time monitoring for reaction cuvettes status, lamp status and temperature. Real time monitoring for reagent volume, reaction curve, calibration curve and QC curve. Linearity range limitation, substrate exhaustion judgement and prozone detecting. Abnormal status warning.	
	Other Function	User authority setting, test panel function, calculated/manual parameters programing, carryover setting, sample and reagent blank auto calculation, automatic failure recovery, automatic print, data statistic, auto/manual dilution test, auto retest.	
	PC Configuration	CPU 2.5GHz, memory 4GB, hard disc 500G 19 inch widescreen monitor	
	Others	PC Interface	RS-232C
Dimension		855(L)x 550(D) x 580(H)	
Power Supply		AC 100-240V, 50/60Hz ± 1 Hz, ≤350VA	
Net Weight		75Kg	
Optional Parts		Barcode reader, ISE module, Water purification module, PC with touch screen, HP laser printer	

AU240-PW

Auto Chemistry Analyzer

Precise, Fast and Stable



 Shenzhen Heto Medical Tech Co., Ltd.

Address: 4/F, Building No.4, Shenzhen Environmental Protection Industrial Park,
No.31 Makan South Road, Nanshan District, Shenzhen, 518055, P.R.China.

Tel: +86 755-28935967

Fax: +86 755-28937587

Email: Sales@hetomedical.com

Web: www.hetomedical.com

 **Contact Us**
Local Distributor

User Friendly Software

- Easy operate with iconic interface
- Multi-language available
- Easy reagent parameters setting
- Support batch request, panel/calculation test
- Print out report template editable
- Auto troubleshoot failure and recovery
- Powerful statistic system
- Manual dilution, auto dilution and post dilution
- Easy rerun and restore test

Calibration & QC

- Multiple calibration mode
- Advanced algorithm to obtain best calibration curve
- Multi QC rules: Westgard, Cumulative sum check, Twin plot



Real Time Monitoring

- Reaction curve, calibration test and QC test curve
- Sample, reagent and reaction cuvette status
- Lamp intensity, water tank and waste tank status
- Reaction disk and reagent disk temperature

AU240-PW

Auto Chemistry Analyzer



Reaction System

- 81 high light transmittance reaction cuvettes
- Metal thermostat ensures stable temperature: 37°C
- Auto water blank test ensures precise results
- Optional quartz glass cuvette



Sample & Reagent Disk

- 24 hours reagent refrigerator: 4-12°C
- Noise free water cooling design
- 40 positions for reagent
- 40 positions for sample with STAT priority
- Optional barcode scanner available



Sampling System

- High polished nano material sample/reagent probe
- High accuracy ceramic syringe
- Auto depth adjustment
- Liquid level detection, vertical and horizontal collision protection
- Maintenance free rotor with stable movement



Anti-contamination Design

- 6-step wash station washes cuvettes with wash solution
- High polished nano material mixer
- Automatic probe and mixer washing
- Intelligent anti-contamination program




Optional

- Internal PC with touch screen
- ISE Module: K^+ , Na^+ , Cl^-
- Barcode Reader

Fluid System

- Branded components including: IWAKI pump, SMC valve, Ceramic syringe, THOMAS peristaltic pump, Tygon tubing and CPC tubing couplings
- Built-in degasser ensures precise sampling

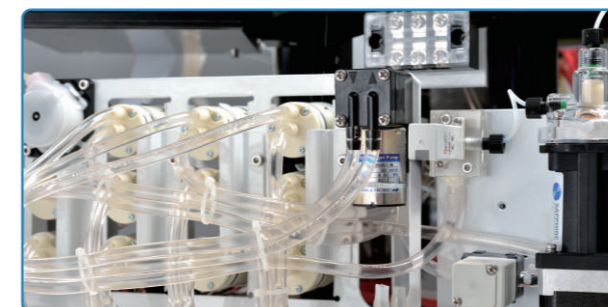

Optical System

- Halogen tungsten lamp with about 2000H life time
- Full-sealed matrix spectrometric system
- Simultaneously dual wavelength detection to avoid interference
- Silent water cooling design

