

Technical Specifications

General Information

Throughput	Up to 480 tests/hour
Measurement Principle	Micron superparamagnetic particles platform With alkaline phosphatase (ALP) labeled reagents And AMPPD substrate

Sample Assembly

Sample Feeder	Continuous sample loading and offloading By sample racks
Sample positions	10 positions per rack, up to 30 racks in one batch
Sample Probe	Steel probe, liquid level detection, clot detection, Horizontal and vertical collision protection
Sample Volume	10~200ul, with increment of 1ul
Sample Dilution	Automatic dilution with 1:2~1:80

Reagent Assembly

Reagent Disk	36 reagent positions, continuous loading, real-time mixing for magnetic particles
Cooling System	2~8 °C
Reagent Pack	50 tests/pack, 100 tests/pack
Reagent Scanning	Built-in/External bar code reader
Reagent Probe	Steel probe, liquid level detection, clot detection, Horizontal and vertical collision protection
Reagent Volume	10~200ul, with increment of 1ul

Substrate System

Substrate Bottle	300/500 tests per bottle
Substrate Loading	Capable of continuous loading during testing
Substrate Heating	Constant heating
Aspiration Volume	200ul

Magnetic separation unit

Separation Assembly	4-phase magnetic separation
Temperature Accuracy	37±0.3 °C

Measurement and Reaction System

Detection Mode	Photon counting
Signal Detector	Photomultiplier (PMT) photo counter
PMT Calibration	LED reference module
Temperature	Accuracy: 37±0.3 °C
Mixing Unit	Non-touch vortexer mixing
Cuvette	Disposable plastic cuvette
Cuvette Loading	Continuous loading and offloading

Working conditions

Dimension	2150mm*1166mm*1300mm (W*D*H)
Weight	580Kg
Noise	<65dBA
Altitude Height	-400m~4000m

CL-6000i

Chemiluminescence Immunoassay system

more than fast




Von-Langen-Weg 10
48565 Steinfurt / Germany
Tel.: 02551/2097
Fax: 02551/80883
E-Mail: info@keul.de
Internet: www.keul.de

www.mindray.com

P/N:ENG-CL-6000i-210285XBP-20191017
©2018 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved.

mindray
healthcare within reach

Fast because of Intelligence

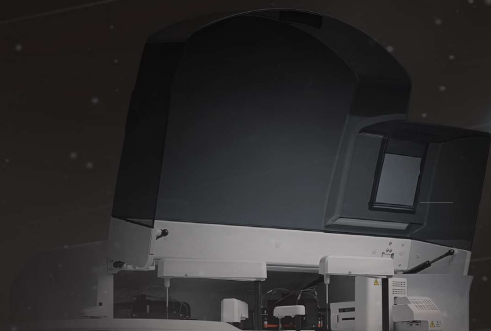




mindray

CL 6000i

Consumables real-time monitoring



All consumables can be loaded
continuously without interruption

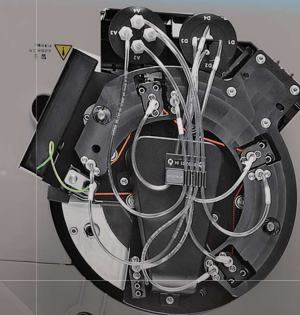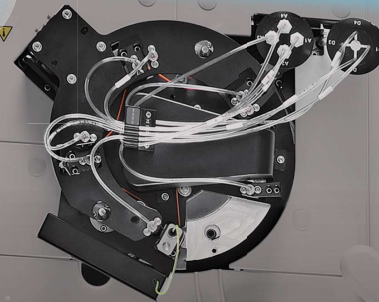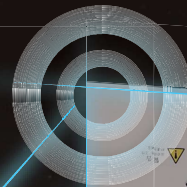

Smart dual-core drive



Doubling the processing power



One instrument serves as two




Up to 480 tests/hour

Fast because of Fluency




Fully automated maintenance

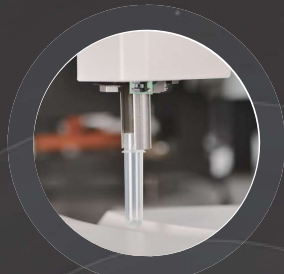

Liquid-solid waste separation
More environmental friendly




Uniquely designed buffer disk to
reduce cuvette jam

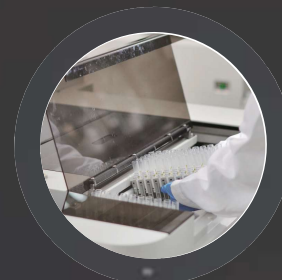

Up to 300 samples can be loaded and
offloaded continuously



2~8 °C refrigerated reagent
carousel featured 36 positions

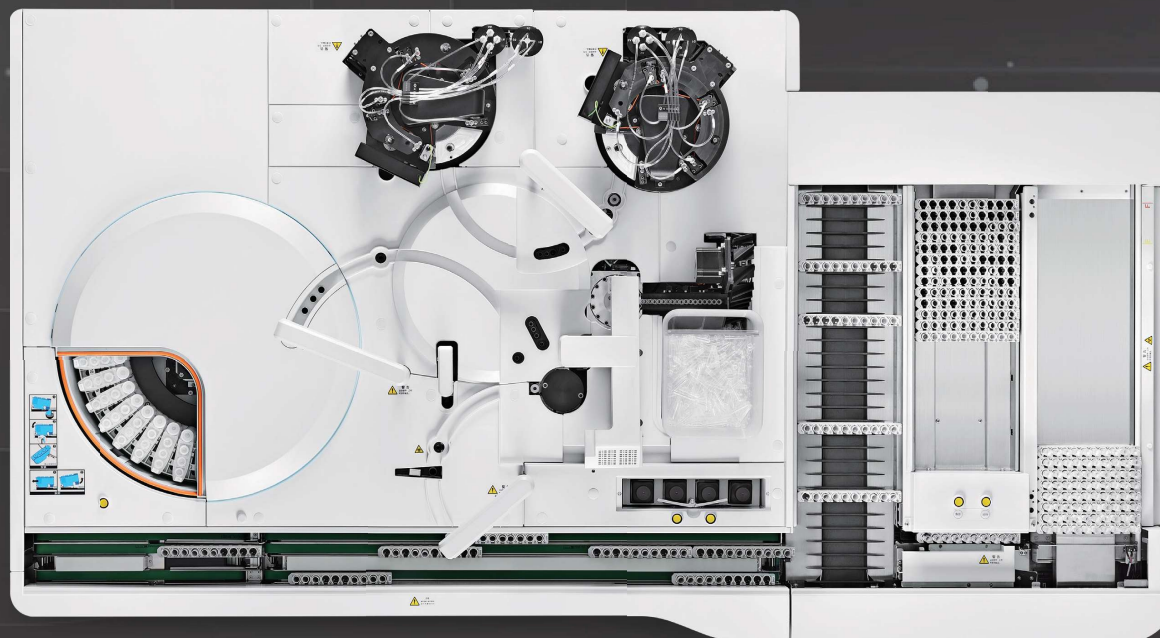

Test Menu

Thyroid

FT3
FT4
T3
T4
TSH
Anti-TG
Anti-TPO
Tg
rT3
TRAB*

Bone Metabolic

Intact PTH
Calcitonin
Vitamin D total

Adrenal Function

DHEA-S
Cortisol
ACTH

Hypertension

Angiotension I*
Angiotension II*
Renin*
Aldosterone*

DS' Screening

PAPP-A*
Free β -HCG*

Infectious Disease

HIV Combo
HBsAg
Anti-HBs
HBeAg
Anti-HBe
Anti-HBc
Anti-HCV
Anti-TP (Syphilis)
HBeAg quantitative*
Anti-HBc IgM*

Cardiac

CK-MB
Myoglobin
Troponin I
BNP

Liver Fibrosis

LN
HA
PIIINP
CIV

Growth Hormone

hGH*
IGF-1*

Fertility

Total β -HCG
FSH
LH
Prolactin
E2
E3
TESTO
PROG
Free testosterone*
17-OH PROG*
SHBG*
AMH*

Diabetes
Insulin
C-peptide

ToRCH

Toxo IgG*
Toxo IgM*
Rubella IgG*
Rubella IgM*
CMV IgG*
CMV IgM*
HSV-1/2 IgG*
HSV-1/2 IgM*
HSV-1 IgG*
HSV-2 IgG*

Tumor marker

CEA
AFP
CA125
CA15-3
CA19-9
FPSA
TPSA
NSE
CYFRA 21-1
CA72-4
PG I
PG II
SCCA
HE4
ProGRP
CA50
CA242

Anemia

Ferritin
Vitamin B12
Folate
RBC Folate

Inflammation

PCT

*: in development

